




## Some ideas how to play with the Flow-Ball KIDS®




Blow the ball up in the air and try to keep it in the air whilst turning yourself round in a circle.



Put a bowl or similar object on the table. Stand back from the table. Take a deep breath, blow through your Flow-Ball KIDS® and try to get the ball into the bowl. Play with a friend and keep score.

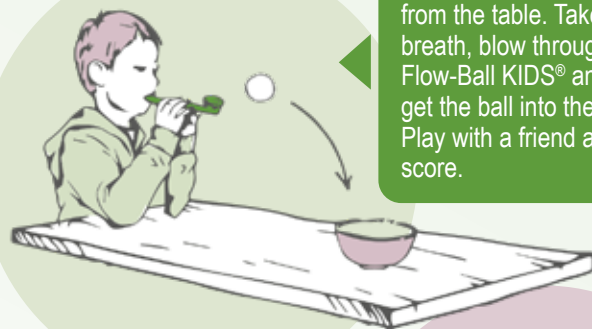

Tilt the Flow-Ball KIDS® slightly upwards or downwards (from horizontal) and blow through it so that the ball spins round at different speeds.

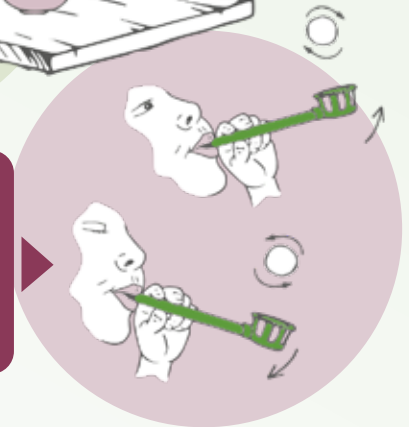


Blow the ball up in the air as high as you can and keep it there for as long as possible.



2 Player Game / Exercise: Tilt the Flow-Ball KIDS® down slightly and blow the ball up high. Your partner should start blowing through their Flow-Ball KIDS® and try to scoop the floating ball from you in mid-air. Now take a breath and try to scoop it back. The aim is to pass the ball backwards and forwards between you without the ball touching either of your baskets or dropping it!



Blow the ball in the air and catch it in the basket, repeat as often as you like.

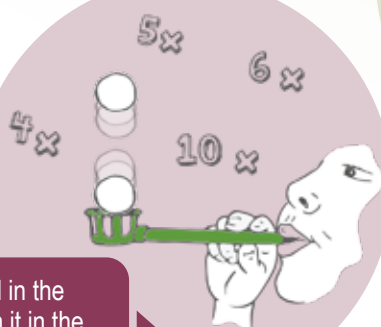

## Care and Maintenance

Because your Flow-Ball KIDS® will be exposed to saliva during use, we recommend that your Flow-Ball KIDS® is cleaned after use. Remove the basket from the stem and wash all parts in soapy water, rinse with clean water, dry thoroughly and reassemble. Once a week we recommend that the device is cleaned using approved cleansing tablets.

The Flow-Ball KIDS® is tested in accordance with EN 71 safety directives.

---

### Head Office:

POWERbreathe International Ltd. (UK REP)  
Northfield Road, Southam,  
Warwickshire, CV47 0FG,  
England, UK.  
Tel: +44 (0)1926 816100  
[www.powerbreathe.com](http://www.powerbreathe.com)

**EC REP** HaB GmbH,  
Porschestra. 4,  
D-21423 Winsen an der Luhe,  
Deutschland.



Flow-Ball KIDS® is not a toy. This product is designed to be used for breathing exercises only. Any other use is not recommended. The material in this manual is for information purposes only. POWERbreathe International Ltd makes no representations or warranties with respect to this manual or with respect to the products described herein. No claims are made or implied in the use or results by the use of equipment herein. POWERbreathe International Ltd shall not be liable for any damages, direct, indirect or incidental, consequential or special, arising out of or related to the misuse or deviations from the intended use of this material or the products described herein.

Flow-Ball KIDS® is a design protected product. All rights reserved. Specifications may change without notice due to manufacturers continuous programme of development.

The POWERbreathe logo type is a registered trademark of POWERbreathe Holdings Ltd. All POWERbreathe product names are trademarks or registered trademarks of POWERbreathe Holdings Ltd.